

Hernik



Rehab Car Seat **STARFIGHT**

... is flexible, comfortable and safe

5-point harness with quick release

reduction inlet: seat depth only 25 cm!

stable thorax supports



Rehab Car Seat **STARFIGHT**

NXT

...flexible, comfortable and safe!

The STARFIGHT-NXT child car seat was developed for children with special needs. It is suitable for installation in most commercial vehicles. The car seat is built for use in the back seat of a car or passenger seat and requires a functioning 3-point safety belt.

In addition to the 3-point belt of your car, the STARFIGHT-NXT is equipped with its own 5-point harness, which provides additional support for your child. This can be adjusted together with the headrest in height.

The back angle is fixed at 90°, however, the car seat can be equipped with a slight (0° - 10°) or extended angular adjustment (10° - 20°), which is often very important for children with limited torso and head control.

The STARFIGHT-NXT is specially designed for very small children, which is why a reinforced two-piece back pad reduces seat depth.

By installing our XL kit, however, the car seat can also be optimally used for later growth - see ②

A novelty in this class are the optional **thorax supports** made of aluminum, which guarantee the child a stable sitting position and are flexibly adjustable.





A very popular feature is the innovative swivel, which facilitates the transfer of the ever heavier child. The pivot was moved from the center to the rear and to the outside. This has been achieved that the seat - during the rotational movement - erases from the backseat of the car and simultaneously moves towards the entrance. Great emphasis was placed on a compact construction: the overall height of the seat is only 16 mm larger!

1

STARFIGHT-NXT grey, with swivel

2

STARFIGHT-NXT, red with optional XL-kit and thorax supports

3

stepless variable and independently adjustable thorax supports

4

ISOFIX adapter for a stable position in the car



...varied, as the needs of their child

The rehab car seat STARFIGHT-NXT has a wide range of optional accessories to meet the child's individual needs. We offer the following:


- additional angular adjustment(10°; 20°)
- swivel (left or right)
- thorax supports made of aluminum
- ISOFIX adapter
- security tray
- abduction bloc
- foot rest
- security pad for the buckle

It is important that no accessory excludes another. Many modern vehicles are equipped with ISOFIX, which is why we offer a suitable adapter system to fix the seat optimally in the car. The 5-point harness can now be comfortably strapped with just one movement, which makes life a whole lot easier.

Since we are constantly developing the STARFIGHT-NXT, we are happy to receive suggestions and ideas!



security tray

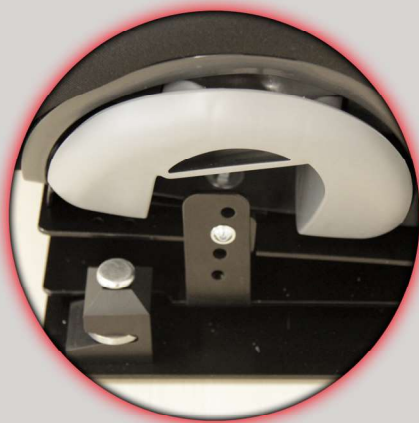

angular adjustment



**Cupholder
to fold-out, on both sides**



**foot rest adjustable
in height and angle**

Your supplier:

HERNIK GmbH
Bodenseestrasse 25
83059 Kolbermoor
☎ +49 (0)8031-2068044
☎ +49 (0)8031-2069215
✉ info@hernik.de
🌐 www.hernik.de